

ORAL HEALTH IMPROVEMENT KIT

to CLEAN & PROTECT YOUR SMILE

GINGIVITIS

DENTISTRX® REVOLUTION® Revolving 360° Toothbrush & UV Sanitizer

- Unique power brush head design effectively reaches around tooth and braces brackets to improve the removal of plaque
- Professionally recommended home-care instrument for disease control
- Built-in timer to guide patients through the ADA recommended two-minute brushing cycle
- UV sanitizer base kills bacteria and viruses on brush head surfaces

Box Set Includes: Brush heads, easy-grip brush handle, power cord, and travel case

[More details at www.DentistRx.com]

PRO-DENRX® STANNOUS FLUORIDE RINSE

- Aids in the prevention of dental caries (cavities)
- Intended for use in the prevention of decay, inhibition of plaque and reduction of tooth sensitivity

[More details at denmat.com]



INSTRUCTIONS FOR ORAL HEALTH IMPROVEMENT KIT

DentistRx®



[Reference: DentistRx® Revolution® owner's manual]

DentistRx® Revolution® requires a tooth-by-tooth brushing technique rather than a scrubbing technique as described here.


- **Brushing Technique:** After applying toothpaste, hold the toothbrush at a 90-degree angle to teeth and gums, and then turn on the power. Allow the toothbrush to clean each tooth by applying light, gentle pressure for a few seconds. Then move to the next tooth and repeat. When using a revolving toothbrush, it is best to follow the gum line and use a “M” pattern when brushing your upper teeth, and a “W” pattern when brushing the lower teeth.
- **Built-In Timer:** The built-in timer will begin the 2-minute cleaning cycle pausing every 30 seconds to indicate when it's time to brush the next section of your mouth.
- **Using the UV Sanitizer:** To sanitize your brush head after use, rinse and shake dry. Place the brush head in the sanitizer tray facing the UV light and close the lid. The sanitizer will automatically run for approximately seven minutes once the lid is closed.

PRO-DENRX® STANNOUS FLUORIDE RINSE

Adults and children 6 years of age and older:

- Use as directed by a dental professional.
- Use once a day after brushing your teeth with a toothpaste or as directed by your dentist.
- Remove cap and seal and replace with pump.
- Rotate spout counter clockwise to extend.
- Depress pump until medicine comes to the tip of the pump.
- Put 1/8 oz. (approx. 2 pump strokes) of concentrated Pro-DenRx Rinse into the provided mixing vial.
- Add water to 1 oz. line on mixing vial
- Tightly snap-on the cap and shake well.

Add water to 1 oz. line on mixing vial

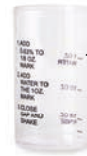

- Vigorously swish between your teeth 1/2 of mixture for one minute and then spit out.
- Repeat procedure with remaining mixture.
- Rinse mixing vial with water after use.
- Do not swallow the rinse.
- Do not eat or drink for 30 minutes after rinsing.

Children 6 to 12 years of age:

- Instruct children in good rinsing habits (to minimize swallowing).
- Supervise children as necessary until capable of using without supervision.

Children under 6 years of age: Consult a dentist or doctor.

[Reference: DenMat product user guide]

Instructions have been provided for patient convenience. Please refer to the manufacturer's product guides included in their package for details.

 **Western Dental®
& Orthodontics**